



Reprocessing Instructions

Kerrison Punches

The following instructions are applicable for Kerrison Punches supplied by Mercian, unless stated otherwise with the packaging of the product.

These instructions conform to the Medical Directive 93/42/EEC (Annex 1 section 13)

These instructions are intended for use only by persons with the required specialist knowledge and training.

Precautions:

- Follow instructions and warnings as issued by manufacturers of any decontaminants, disinfectants and cleaning agents used. Wherever possible avoid use of mineral acids and harsh, abrasive agents.
- **Do not** clean instruments in Theatres with saline, use sterile water.
- Visually inspect the Kerrison, in particular the **location screws to ensure they are secure.**
- The shaft & tip will require particular attention during cleaning.
- As per HTM 0101 1.138 alkaline detergents in the pH range **8.0–11.0** are preferred.


Inspection;

- Inspect all instruments prior to sterilisation or storage to ensure instruments are suitable for use. Set aside any instruments that are fractured, bent or have damaged tips; send these instruments for service or repair.

Note: **Do not remove traverse screws during reprocessing.** Doing so can lead to screws backing-out and subsequent instrument failure.



Ensure screws are secured and seated flush.

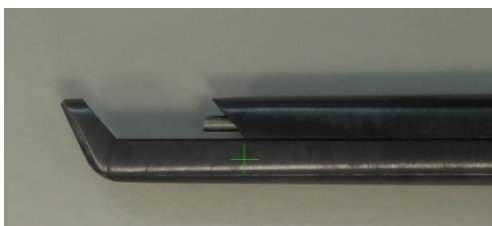

Check tips for damage

Limitations on Reprocessing:

- Repeated processing has minimal effect on these devices.
- End of life is normally determined by wear and damage in use.
- Any specific limitations on the number of reprocessing cycles are indicated in the IFU for the device.

INSTRUCTIONS

Kerrison Punches should only be used for their intended purpose; this Kerrison Punch is used in Lumbar and Cervical Spinal Surgery.

From Point of Use:

- Wherever possible, do not allow blood, debris or bodily fluids to dry on devices. To prolong the life of the device, reprocess as soon as possible after use and within 6 hours of use. If they cannot be reprocessed immediately, use an enzymatic foam spray cleaner to help prevent soil from drying.

Cleaning - Automated washer-disinfector machines:

Our Instruments are suitable to be processed through automated washer-disinfector machines which are CE marked and validated to:

- ISO 15883-2:2009 *Requirements and tests for washer-disinfectors employing thermal disinfection for surgical instruments.*
 - Adhering to: Health Technical Memorandum: Management and decontamination of surgical instruments: Part D - Washer-disinfectors
1. Load devices carefully, with any box joints and hinges open and so that any fenestrations in devices can drain.
 2. Place heavy devices with care in the bottom of containers, taking care not to overload wash baskets.
 3. Place devices with concave surfaces (e.g., curettes) facing down to prevent pooling of water.
 4. Where available, use appropriate attachments to flush inside reamers and devices with lumens or cannula.
 5. Ensure that soft, high purity water which is controlled for bacterial endotoxins is used in the final rinse stage.

Cleaning parameters:

Process:	Parameters:
Washer Disinfector (Thermal)	90°C (-0°C +5°C) for a minimum of 1 minute
Autoclaving	Please refer to the 'Sterilisation' section of these instructions.
Final rise	Tap water / Reverse osmosis water is acceptable

Note: *Automated cleaning may not be suitable for all lumens and cannula, in which case clean manually with a water jet gun, if available, and an appropriate brush (and stilette if provided) that reaches the depth of the feature. After manually cleaning, pass all devices through an automated cleaning cycle to achieve disinfection.*

Cleaning - Manual:

- **Manual cleaning is not advised if an automatic washer-disinfector is available.**
If this equipment is unavailable, use the following process:
1. Using a double sink system (wash / rinse) dedicated for device cleaning (not used for hand washing). Ensure that the water temperature does not exceed 35°C.
 2. In the first sink, keeping the device submerged, with an autoclavable brush, apply CE marked cleaning solution to all surfaces until all soil has been removed. Pay particular attention to serrations, teeth, ratchets and hinges, always brushing away from the body and avoiding splashing. Ensure rongeurs and hinged devices are thoroughly cleaned in both open and closed positions.
 3. In the second sink, rinse the device thoroughly with soft, high purity water which is controlled for bacterial endotoxins, so that the water reaches all parts of the device, then carefully hand dry or use a drying cabinet.

Note: *Manual cleaning is NOT a disinfection process. When manual cleaning is used it may not be possible to disinfect the device prior to further handling.*

Cleaning - Inspection:

- After cleaning, visually inspect all surfaces, cannulations, ratchets, joints for complete removal of soil and fluids. If any soil or fluid is still visible, return the device for repeat decontamination.

Maintenance:

- Apply surgical grade lubricants to hinges, joints and moving parts as per the lubricant manufacturer's instructions.

Inspection & Function Testing:

- Visually inspect and check all devices for damage.

Note: *If a device is returned to the manufacturer/supplier, the device MUST be decontaminated and sterilised and be accompanied by the relevant documented evidence.*

Packaging:

- All devices are to be packed following local protocol in accordance with BS standards.

Sterilisation:

- Sterilize in a steam autoclave conforming to BS EN 285:2015 at a holding temperature of 134°C to 137°C for between 3 to 3.5 minutes.
 - Validated to the following standard EN 17665-1:2006
 - Adhering to: Health Technical Memorandum 01-01: *Management and decontamination of surgical instruments*. Part C: Steam sterilization

Storage:

- Ensure devices are dry before storage, and stored in dry, clean conditions at an ambient room temperature.

Additional Information:

- Other forms of cleaning (i.e., ultrasonic) and sterilisation (i.e., low temperature steam and formaldehyde, ethylene oxide and gas plasma) are available. However, always follow the instructions for use issued by the manufacturer and always consult with them if in any doubt over the suitability of any process used.
- Cleaning and sterilizing guidelines are detailed in:
Health Technical Memorandum 01-01: *Management and decontamination of surgical instruments*.
Part C: Steam sterilization Part D: Washer-disinfectors.
<https://www.gov.uk/government/publications/management-and-decontamination-of-surgical-instruments-used-in-acute-care>

Note: *It is the responsibility of the reprocessor to ensure that the reprocessing as actually performed using equipment, materials and personnel in the reprocessing facility achieve the desired result. This requires validation and routine monitoring of the process. Likewise, any deviation by the reprocessor from the instructions provided must be properly evaluated for effectiveness and potential adverse consequences.*



Mercian Surgical Supply Co Limited
10 Topaz Business Park
Topaz Way
Bromsgrove
UK
Tele: 0844 879 1133
Email: info@merciansurgical.com