


# The miniPrevail<sup>®</sup> Family

Your artistry goes hand in hand with our design.

For the masterpiece you create in every smile.

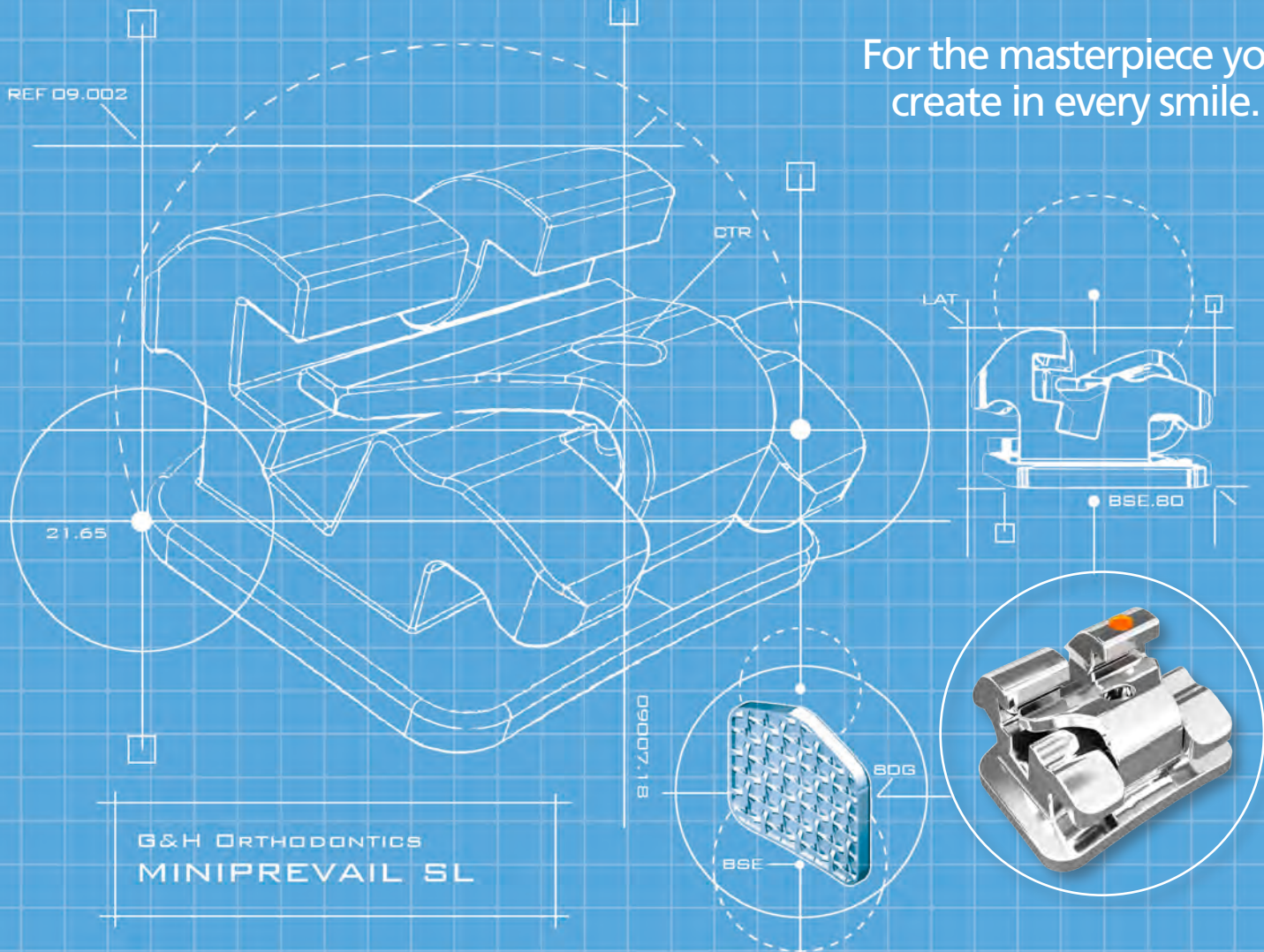

## miniPrevail<sup>®</sup>TWIN

TOP SELLER



A low profile, miniature twin bracket that incorporates all the benefits of a larger bracket.

## miniPrevail<sup>®</sup>TUBE



Our distinct design takes the guesswork out of placement while the tube features offer easy wire placement and engagement.

## miniPrevail<sup>®</sup>SL



An interactive self-ligating bracket. The answer to achieve excellent, predictable results in less time.

Shop our full line of G&H manufactured products @ [GHOrthodontics.com](http://GHOrthodontics.com)

• U.S. & Canada Toll-Free: 800-526-1026 • International +1 317-346-6655



ORTHODONTICS<sup>®</sup>  
Quality Manufacturing. Exceptional Service.



TOP SELLER

# miniPrevail<sup>®</sup> TWIN

LOW PROFILE MINIATURE TWIN BRACKET SYSTEM

## Our low profile bracket precision manufactured by G&H<sup>®</sup> Orthodontics

miniPrevail<sup>®</sup> TWIN brackets are constructed of a specialized micro-grain stainless steel MIM powder, ideal for creating small, robust appliances that are inherently stronger and more dimensionally accurate than other similar brackets.

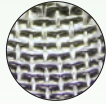Intelligent Design	Patient Comfort
<p><b>METAL INJECTION MOLDED (MIM)</b> micro-grain stainless steel powder ensures exact tolerances and strength of design</p> <p><b>SURGICAL-GRADE STAINLESS STEEL</b> for greater strength</p> <p><b>COLOR ENHANCED I.D.</b> for easy bracket selection</p> <p><b>CENTER SCRIBE LINE</b> aligns with long axis for accurate placement</p> <p><b>TORQUE IN BASE</b> for level slot lineup and improved finishing</p> <p><b>80 GAUGE ANATOMICAL CONTOUR BONDING BASE</b> matches the anatomical curvature of the tooth for maximum contact and a strong, consistent bond</p> 	<p><b>LOW PROFILE</b> • reduces occlusal interference • enhances comfort</p> <p><b>ROUNDED FACIAL CONTOURS</b> for patient comfort</p> <p><b>INCREASED TIE WING AREA</b> to facilitate ligation</p> <p><b>RHOMBOIDAL SHAPE</b> Mesiodistal sides line up with the long axis and incisogingival sides line up with the occlusal plane</p> 

### Did you know?

Our **ELASTOMERICS** are engineered to exacting standards to ensure consistent force delivery. G&H elastomerics, archwires, and tubes are designed to work seamlessly with our bracket systems.

## miniPrevail<sup>®</sup> TWIN

Available in .018 or .022 slot size  
Roth\* or MBT\* prescription

Sold individually or in kits.

1 kit = 20 brackets oriented in a kit box.



Also available in Roth\* V-Slot  
(Vertical slot measures .018" x .018")  
Sold individually or in kits.

Shop our full line of G&H manufactured products @ [GHOrthodontics.com](http://GHOrthodontics.com)

# miniPrevail<sup>®</sup> SL

INTERACTIVE SELF-LIGATING BRACKET SYSTEM

**The answer to achieve excellent, predictable treatment results in less time.**

An **INTERACTIVE** self ligating bracket which uses a NiTi spring clip that is strong enough to deliver continuous and consistent forces while assuring tooth movement control. When using round archwires (up to .020") allows **PASSIVE** engagement during initial leveling and alignment to reduce friction and improve sliding mechanics. When using rectangular archwires (from .016"x.022"), becomes **ACTIVE** to provide effective torque and rotation control toward the end of treatment.

## Intelligent Design

**COLOR ENHANCED I.D.**  
for easy bracket selection

**ASSEMBLED USING LASER TECHNOLOGY**

**17-4 SURGICAL-GRADE STAINLESS STEEL**  
for greater strength

**METAL INJECTION MOLDED (MIM)**  
Constructed of a specialized micro-grain stainless steel MIM powder, ideal for creating small, robust appliances that are inherently stronger and more dimensionally accurate than other similar brackets. Ensures exact tolerances and strength of design

**COMPOUND CONTOUR BASE**  
for increased bond strength and precise bracket placement



**MICRO-ETCHED 80 GAUGE MESH ANATOMICAL BONDING BASE**  
matches the anatomical curvature of the tooth for maximum contact and a strong, consistent bond

Also –

Take a look at our **PREMIUM ARCHWIRES**, we make our wires, and that makes the difference. Premium alloys, material science and engineering, precision manufactured in the U.S.A.

See our website or catalog for full details.

## Patient Comfort

**LOW PROFILE**


- reduces occlusal interference
- enhances comfort

**AMPLE TIE WING**

for elastomeric chain and ties

**FULL SLOT AREA COVERAGE**

for total tooth movement control

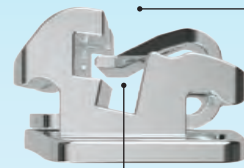

Optional hook on cuspid and bicuspid



**NITI SPRING CLIP**

for maximum resistance to permanent deformation, durability, and stable performance - prevents unwanted opening, yet opens easily with an explorer, closes with gentle finger pressure

## miniPrevail<sup>®</sup> SL

Available in .018 or .022 slot size  
Roth\* or MBT\* prescription

Sold in packages of 5 or in kits.

1 kit = 20 brackets oriented in a kit box.



Visit [GHOrthodontics.com/miniPrevailSL](http://GHOrthodontics.com/miniPrevailSL)

**Not for sale in USA or Canada**

Shop our full line of G&H manufactured products @ [GHOrthodontics.com](http://GHOrthodontics.com)

# miniPrevail<sup>®</sup>TUBE

PRECISION MANUFACTURED BY G&H<sup>®</sup> ORTHODONTICS

## Seamlessly working together with our miniPrevail<sup>®</sup>TWIN and SL

Our new miniPrevail<sup>®</sup>TUBE has been engineered with the patient and orthodontist in mind. Our distinct "crosshair" design takes the guesswork out of placement while the tube features offer easy wire placement and engagement. Finally, the low profile and sleek hook is designed for the utmost patient comfort.

Intelligent Design	Patient Comfort
<p><b>DISTINCT "CROSSHAIR"</b> horizontal and vertical placement sight lines with color-coded ID center "bullseye" positioning feature</p> <p><b>OCCUSAL/GINGIVAL PLACEMENT GUIDES</b> allow for easy placement at the posterior region and assist with minor adjustments prior to curing the adhesive</p> <p><b>FUNNELED ENTRANCE</b> for easier wire insertion</p> <p><b>ANATOMICALLY FORMED BASE MICRO-ETCHED 80 GAUGE MESH</b> with developmental indent for precise fit and positioning, micro-etched mesh provides enhanced bond strength</p> <p><b>BALANCED</b> to facilitate easier bonding</p> <p><b>PARALLEL WALLS</b> for improved grip when handling</p>	<p>Available in 4.2mm or 3.2mm <b>MESIAL/DISTAL DIMENSION</b> for full rotational control</p> <p><b>LOW PROFILE</b> metal injection molded (MIM) stainless steel</p> <p><b>ANGLED OCCUSAL WALL</b> to avoid interference</p> <p><b>SLEEK HOOK DESIGN</b> for ultimate patient comfort and simplified bonding</p> <p><b>Mesh Base Pad 1st Molar Tube</b> Width: .260" (6.604mm) Height: .120" (3.048mm)</p> <p><b>Mesh Base Pad 2nd Molar Tube</b> Width: .200" (5.08mm) Height: .1150" (2.921mm)</p>

**Need BANDS?**  
Our **BANDS**, provide a secure fit you can count on. Precision-shaped tooth anatomy for a tight, secure fit and excellent retention, time after time.  
See our website or catalog for full details.

## miniPrevail<sup>®</sup>TUBE

Available in .018 or .022  
Weldable or Bondable  
Roth\*, MBT\* or Standard prescription  
Sold in packages of 10.

Shop our full line of G&H manufactured products @ [GHOrthodontics.com](http://GHOrthodontics.com)

• U.S. & Canada Toll-Free: 800-526-1026 • International +1 317-346-6655

