


# Progeny ResQ

**Reduce fluorescence interference  
with the 1064nm advantage –  
the new generation in  
handheld chemical detection.**

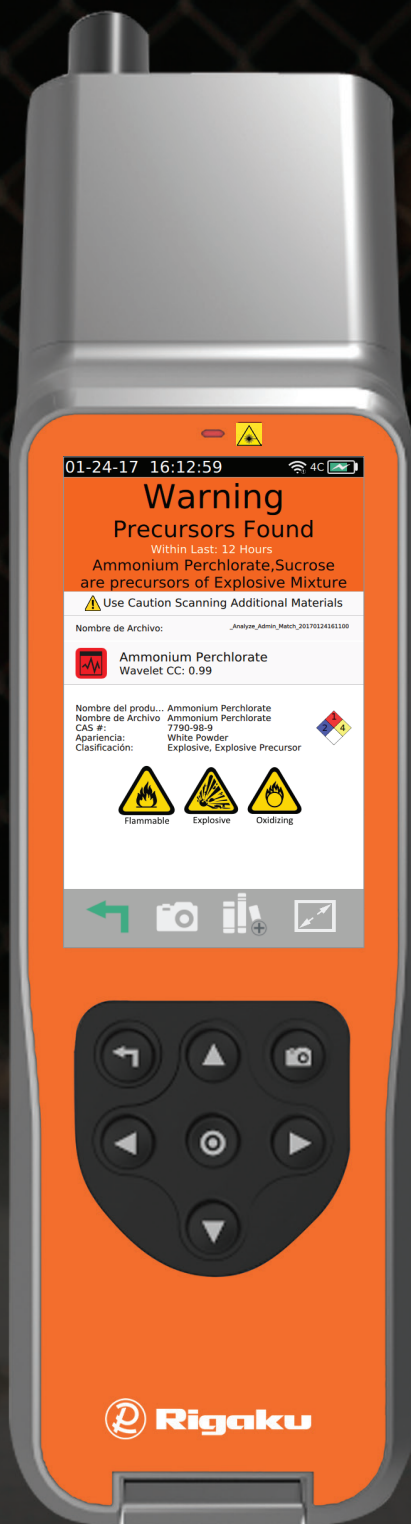

For first responders and law enforcement agencies who are facing an increase in sophisticated threats and the circulation of illegal narcotics, Progeny® ResQ® provides the industry's most comprehensive chemical detection range in a handheld form. With the use of 1064nm Raman, Progeny ResQ reduces sample fluorescence interference issues, allowing emergency responders to identify more substances - including colored materials or through containers.



By integrating Progeny ResQ into routine response tactics, users can expect the ability to:

- **Collect analysis in any environment with a rugged, MIL-STD-810G tested tool**
- **Easily navigate using large buttons and a bright screen display for use in the field**
- **Increase confidence with largest standard library or easily adapt for emerging threats**
- **Expand incident response with automatic alerts of potentially dangerous combinations**


## PROGENY RESQ SPECIFICATIONS

Mode of Operation	Handheld, with optional docking station	Efficiently identify threats at the point of need Includes measurement delay timer (up to 10 minutes) for additional safety
Excitation Wavelength	1064nm	Measure a broader range of substances
Output Power	up to 490 mW adjustable laser excitation power	Allows for customization of analysis parameters to obtain most accurate analysis of your specific materials
Graphical User Interface	Touchscreen and softkeys	Smartphone-inspired with bright, touchscreen display or large button navigation for protective glove use
Digital Camera	On-board	Record sample evidence and store with analysis
Reporting Capabilities	On-board	Immediate PDF generation and transfer of results with camera images
Included Standard Library	Yes	Up to 13,000 compounds including CWA, explosives, TIC/TIM, narcotics
User-Created Library	Yes	Easily adaptable for new, emerging threats into custom library; export to share with other Progeny ResQ users worldwide
Threat Recipe Alerts with 4C™ Technology	Automatic	Expands incident response by warning users of potential dangerous combinations
Data Transfer	Yes	Transfer results via WiFi or USB
Operating Temperature	-20 to 50°C intermittent, +40°C continuous	
Certification	MIL-810-G, IP-68	
Battery	14.4V Li-ion, 3100mAh, > 5hrs	Long-lasting charge with options for a recharger or docking station for 24/7 operation
External Power	100~240VAC, 50-60Hz, +18V DC, 2.5A	
Size	29.9cm x 8.1cm x 7.4cm (11.8in x 3.2in x 2.9in)	Allows for single-hand operation
Weight	~1.6kg (3.7lbs)	
Warranty	24 months	
Industry Regulation Compliance	SWGDRUG Category A technique for presumptive testing FDA 21 CFR Part 11 compliance	

Progeny ResQ is not ITAR-restricted, though a U.S. Department of Commerce export license may be required for some countries.

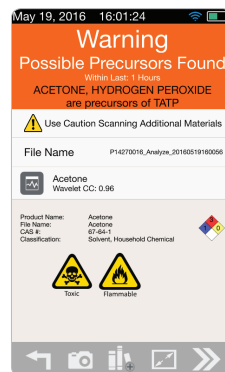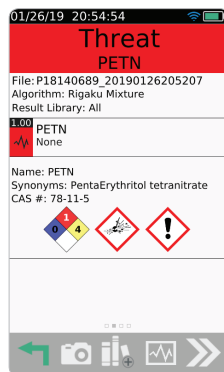

All products manufactured by Rigaku Analytical Devices Inc. are made in the USA. ©2020 Rigaku Analytical Devices, Inc. Progeny and ResQ are registered trademarks and 4C is a commercial trademark of Rigaku Analytical Devices, Inc. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representatives for details.

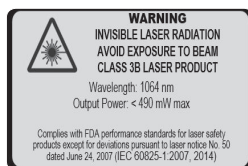

**Rigaku Analytical Devices, Inc.**  
Boston, MA

Toll Free: +1 855.785.1064

Direct: +1 781.328.1024

Email: [handhelds@rigaku.com](mailto:handhelds@rigaku.com)

[www.rigaku.com/ResQ](http://www.rigaku.com/ResQ)